
DK .....	Sejl med multistøtte.....	2
GB/US .....	Multi support sling.....	15
DE.....	Multi-Stützsitz.....	29
SE.....	Multi support sele .....	43
FR.....	Harnais multi support.....	56
IT .....	Imbragatura multi support.....	70

User manual – vers. 101.00

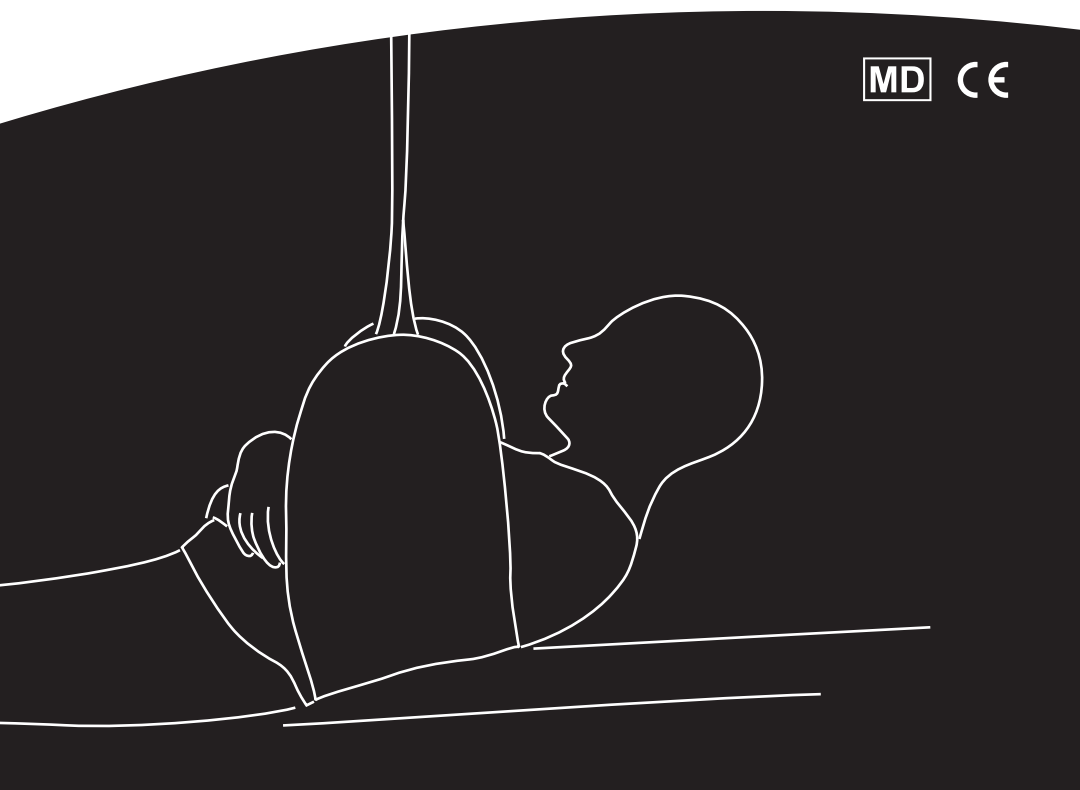Guldmann™

# DK . . . . . Sejl med multistøtte

Vers. 101.00

Varenr.:

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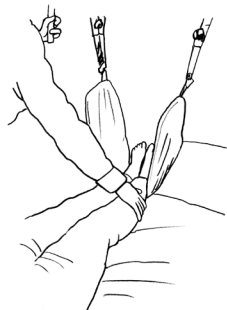
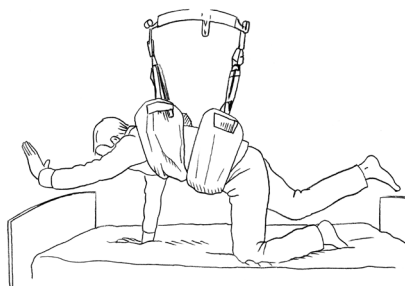
Vers. 101.00

Item nos:

28467 Multi support sling, yellow polyester  
284223 Multi support sling, white polyester

---

<b>1.00</b>	<b>Purpose and use</b>	<b>16</b>
1.01	Manufacturer	16
1.02	Intended purpose	16
1.03	Area of use	16
1.04	Conditions of use	16
1.05	Important/Precautions	16
1.06	UHF RFID tag	17
1.07	Labels and Marking	17
1.08	Use	18
<hr/>		
<b>2.00</b>	<b>Maintenance</b>	<b>19</b>
2.01	Cleaning	19
2.02	The owner's daily maintenance duty	19
2.03	Disposal of slings	20
<hr/>		
<b>3.00</b>	<b>Service and lifetime</b>	<b>20</b>
3.01	Safety/service inspections	20
3.02	Lifetime	22
<hr/>		
<b>4.00</b>	<b>Technical specifications</b>	<b>23</b>
<hr/>		
<b>5.00</b>	<b>EU-Declaration of conformity</b>	<b>23</b>
<hr/>		
<b>6.00</b>	<b>Environmental policy statement – V. Goldmann A/S</b>	<b>23</b>
<hr/>		
<b>7.00</b>	<b>Warranty and service conditions</b>	<b>23</b>
A.	Warranty	23
B.	Service or Repair	24
<hr/>		
<b>8.00</b>	<b>Fitting Multi support sling – lifting upper body</b>	<b>25</b>
<hr/>		
<b>9.00</b>	<b>Fitting Multi support sling – lifting leg from bed</b>	<b>26</b>
<hr/>		
<b>10.00</b>	<b>Fitting Multi support sling – lifting leg from wheelchair</b>	<b>27</b>
<hr/>		
<b>11.00</b>	<b>Fitting Multi support sling – other uses</b>	<b>27</b>

---

**1.00 Purpose and use**

---

**1.01 Manufacturer**

V. Guldmann A/S  
Graham Bells Vej 21-23A  
DK-8200 Aarhus N  
Tel. +45 8741 3100  
www.guldmann.com

---

**1.02 Intended purpose**

The sling is intended for lifting or supporting a person or body parts of a person.

---

**1.03 Area of use**

The sling is suited for use in hospitals, nursing homes, institutions, rehabilitation centers and in private homes.

---

**1.04 Conditions of use**

The sling is designed for use in both mobile lifters and lifting hoist systems. The multi support sling is ideal for lifting various body parts, such as heavy legs in connection with fitting lifting slings, or when personal hygiene and other patient care tasks are performed.

The sling can be used to lift the patient's upper body, lower leg, foot, arm and pelvis.

The following prerequisites apply whenever the sling is used:

- The sling must only be used by qualified personnel or individuals who have received instruction in how to use the sling in question.
  - The correct sling size must be used.
  - The maximum lifting capacity of 225 kg must never be exceeded.
  - The caregiver must be aware of the patient's well-being when using the lifting sling.
  - The sling is to be used with a Guldmann lifting hanger.
- 

**1.05 Important/Precautions**

- Read the instructions carefully before using the sling.
- The slings maximum load must never be exceeded.
- The sling may only be used to lift a person.
- Before a sling is used, it must be examined according to point 2.02.
- Never use a sling that is too big for the user.
- Possible repairs must only be made by the manufacturer.
- Any serious incident that occurred in relation to this device should be reported to the manufacturer and the local competent authority.

1.06

UHF RFID tag

This product features a passive UHF RFID tag integrated in the product label. The RFID tag can be used for asset management & tracking purposes. The RFID tag is readable with equipment compliant to EPC global UHF Class I Gen ISO 18000-63.

1.07

Labels and Marking

CE marking


Medical Device Class I in accordance with EU MDR Regulation


Read the manual before use


Example of serial number label

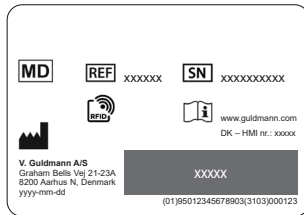
Inspection label

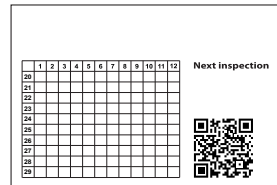Product label

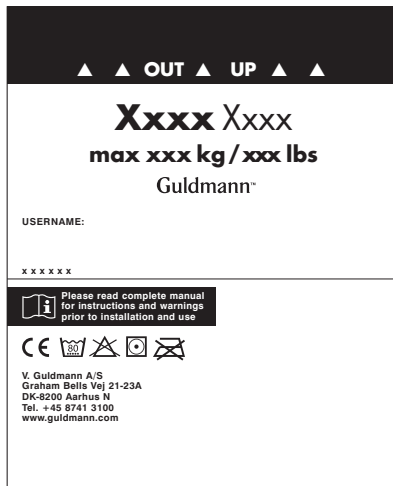Size label

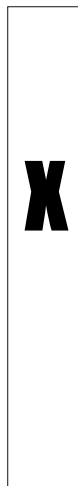**Use**

***If there is any doubt about the selection or use of a lifting sling, please contact your supplier.***


**Important!**

Plan the moving procedure and do not use the sling without supervision. Before lifting, ensure that the patient cannot become stuck and that the sling is not caught on the bed, wheelchair or on other surrounding items. The patient's head, arms, hands and feet must not be in danger of becoming trapped. Be careful with any tubes or wires that may be attached to the patient and/or nearby equipment. Check that the hand control and hand control cable are free of the lifting hanger, patient and other objects before the lifting hanger is moved up or down.


Guldmann shall not be liable for faults or accidents due to incorrect use of the lifting sling, or for reasons of inadequate attention on the part of the carer or user. If the sling is used in combination with products that are not manufactured by Guldmann, a risk assessment must be made by qualified staff.

**Lifting hanger, 4 attachment points****Caution!**

**Be careful when attaching the lifting sling's straps on the hooks. Check that the straps have been correctly placed in the lifting hanger's hooks. When pressing the up button on the hand control to lift the user, check again that all straps remain correctly placed in the lifting hanger's hooks (Fig. 1).**

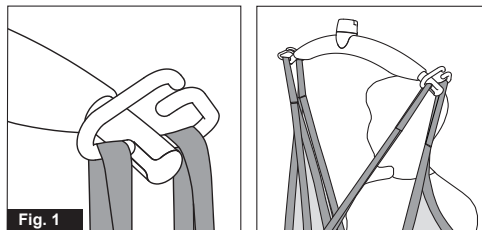**Lifting hanger****Caution!**

**Be careful when attaching the lifting sling on the hooks. Check that the straps have been pulled completely through the rubber safety catch (A) and into place in the lifting hanger's hooks. When pressing the up button to lift the user, check again that all the straps remain correctly placed in the lifting hanger's hooks (fig. 1a and fig. 1b).**

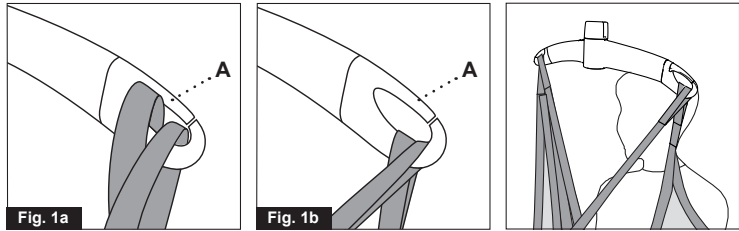Placing the sling, look at section 8

---

## 2.00 Maintenance

### 2.01 Cleaning

Normal washing at the indicated temperature

**Note:** The sling in white material can be washed max. 5 times

Do not use bleaching agent

Tumble-drying at low temperature

Do not iron

---

### 2.02 The owner's daily maintenance duty

Check the lifting sling for wear and damage before use according to the following checklist which is not intended to represent all potential inspection steps. Potential damage may vary. Judgment of inspector/site prevails.

#### Sling inspection checklist

Before using a Guldmann sling / accessory check the following:

##### Is the sling clean?

Follow facility specific infection control procedure.

##### Is the sling's label present, legible and complete?

Missing, illegible or incomplete sling label(s) could make identification of appropriate size of the sling, function of sling, and or weight limit capacity of the sling impossible.

##### Are the lifting straps and stitches intact?

- Look for broken or worn stitches
- Look for knots in straps
- Look for tears or fraying of straps
- Look for snags or punctures or holes
- Look for any particles in fabric or straps

**Is the fabric intact?**

- Look for abnormal wear patterns, excessive wear, abrasive evidence
- Look for cuts or frayed fabric
- Look for unusual or significant discoloration
- Look for snags, punctures, tears, holes
- Look for frayed or insecure seams
- Look for any acid / caustic / thermal burns
- Look for changes in material consistency, e.g. increased stiffness
- Look for any imbedded particles

**Has the shape of the sling been altered, made shorter or longer in relation to the original size using knots, needles, tape or other methods?**

**Conclusion**

If the sling suffers from one or more of the above mentioned conditions then it must be taken out of service regardless of the weight of the person to be lifted.

---

**2.03**

**Disposal of slings**

Slings are disposed of by incineration. By proper incineration polyester will be degraded to carbon dioxide and water.

---

**3.00**

**Service and lifetime**

---

**3.01**

**Safety/service inspections**

In accordance with international standard EN/ISO 10535 “Hoist for the transfer of disabled persons – Requirements and test methods” an inspection **must** be performed every 6-month according to the following instructions, which is not intended to represent all potential inspection steps. Potential damage may vary. Judgment of inspector/site prevails.

**Safe Operating Practices with Slings**

Considerations for damaged or defective slings and taking them out of service:

Withdraw the sling from service if one or more of the following conditions are present:

- chemical or caustic burns
- melting or charring of any part of the sling
- snags, punctures, tears or cuts
- broken or worn stitches
- missing, illegible or incomplete sling tag
- knots in any part of the sling
- abrasion
- other visible damage that causes doubt as to the strength of the sling

Sling inspection is done for the protection of the user, the caregiver, and the overall hospital site safety. A sling inspection system has additional benefit. Systematic sling inspection will assist in the identification of damage trends, potentially leading to cost effective suggestions and results. The inspection process can also help to identify inventory duplicity in certain sling types and sizes.

### **Sling inspection system**

Development of a specific procedure and program for the inspection of slings at your facility is your best safeguard. Consider employing a three part system of inspection. Slings that are removed from service and are not capable of repair should be disposed of so they are unfit for any future use and can not find a way back into active inventory.

#### **1) Initial**

This level of inspection is done at the time that the sling is received into your facility. The inspector should ensure that no damage has occurred during transit, and also verify that the sling work load limits match those contained in the manufacturer's catalogue. If your facility documents the sling inspection process through written inspection records, the paper trail should begin at this stage

#### **2) Frequent**

The frequent level of inspection should be done by the sling user before each use. The sling should be examined and removed from service if damage is detected. The sling user should also determine that the sling is proper for the user conditions, care task required and the required weight capacity.

#### **3) Periodic**

Your facility might want to consider implementing a program for a periodic level of inspection at regular intervals. The interval should be based upon the frequency of use, severity of the service cycle and information derived through the inspection process. Recommendations to prevent damage and enhance service life could be made by staff that perform the periodic inspections. If written inspection records are maintained, they should always reference the unique sling identification number, and be updated to record the condition of the sling. Not intended to represent all potential inspection steps or all potential aspects of product management program. Judgment of inspector/site prevails.

### **Sling inspection technique**

The sling inspection procedure should be thorough, systematic and consistent; both visual and "hands on" inspection techniques are recommended. Certain forms of damage are far more discernable through hands-on inspection, than by visual inspection. For example, fabric stiffness, crushed webbing, as well as, thinning fabric can be identified through tactile inspection. Visual inspection alone may not reveal all forms of sling damage. Once signs of damage have been identified, do not downgrade the work load limit of the sling, with the intent of continuing to use it, but at limited capacity or frequency. This is sometimes done to get more service life out of a damaged sling. The operating rule and standard should be: intact = use; damage = do not use.

Consider the practice of documenting sling inspections through written inspection records. The documentation should include information such as: the name of manufacturer, the sling stock number, width and length, the unique sling identification number (important in differentiating similar slings), as well as the condition of the sling. Other important information might also include the date it was received or put into use at your facility and any special features (if applicable). A beneficial outcome of an inspection program would be the realization of repetitive forms of damage and the analysis that would lead to specific recommendations.

**Sample visual examples of synthetic sling damage <sup>x)</sup>**


*Chemical/caustic burns*

*Broken stitching*

*Crushed / Frayed webbing*

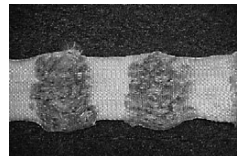*Knots*

*Melting / Charring*

<sup>x)</sup> sample visual images not intended to represent all types of potential damage

**3.02**

**Lifetime**

The life of the sling is individual and depends on how it is used, washed etc. Before use the sling must be examined according to description in section 2.02 and if it does not meet the inspection requirements, it must be discarded if necessary.

The sling in white material can be washed max. 5 times.  
Remember to indicate on the label every time the sling is washed.  
Discard the sling after 5 washing.

---

**4.00 Technical specifications**

Lifting capacity, SWL . . . . .	255 kg (560 lbs)
Material, 28467 . . . . .	Polyester (yellow)
Material, 284223 . . . . .	Polyester (white)

---

**5.00 EU-Declaration of conformity**

The product is manufactured in compliance with regulation (EU) 2017/745 of the European parliament and of the Council of 5 April 2017, as medical device class I.

---

**6.00 Environmental policy statement – V. Guldmann A/S**

Guldmann is continuously working towards ensuring that the company’s impact on the environment, locally and globally, is reduced to a minimum.

It is Guldmann’s goal to:

- Comply with the current environmental legislation (e.g. WEEE and REACH directives)
- Ensure that we, at the widest possible range, use RoHS compliant materials and components
- Ensure that our products do not have an unnecessary negative impact on the environment regarding use, recirculation or disposal
- Ensure that our products contribute to a positive working environment in the places they are utilised

Inspections are made annually by the Department for Nature and Environment from the Municipality of Aarhus using the Danish Environmental Protection Act, section 42 as a reference.

---

**7.00 Warranty and service conditions**

---

**A. Warranty**

Guldmann warrants its equipment is free from material defects under normal use, and will perform substantially in accordance with the specifications set forth in documentation provided with the equipment.

This express warranty shall be in effect for one year from the date of original purchase and installation (the “Warranty Period”). If a valid claim is made during the Warranty Period for malfunction or equipment defect, Guldmann will repair or replace the equipment at no additional cost to you. Guldmann retains sole discretion as to whether the equipment will be repaired or replaced.

The warranty does not cover any part of the equipment that has been subject to damage or abuse by the user or others. The warranty does not cover any part of the equipment that has been altered or changed in any way by the

user or others. Guldmann does not warrant that the lifting device functions will meet your requirements, be uninterrupted or error free.

The warranty set forth is in lieu of all other express and implied warranties, whether oral, written or implied, and the remedies set forth above are your sole and exclusive remedies. Only an authorized officer of Guldmann may make modifications to this warranty, or additional warranties binding on Guldmann. Accordingly, additional statements such as advertising or presentations, whether oral or written, do not constitute warranties by Guldmann.

This warranty shall be null and void if the equipment is operated and maintained in any manner inconsistent with its intended use or the instructions provided with the product. Further, in order for the warranty to remain in effect for the full Warranty Period, all service to the equipment must be provided by a Guldmann certified technician. Any parts or components repaired or replaced by a Guldmann certified technician will be guaranteed for the remainder of the Warranty Period.

#### **Only for USA**

This warranty shall be null and void if the equipment is operated and maintained in any manner inconsistent with its intended use or the instructions provided with the product. Further, in order for the warranty to remain in effect for the full Warranty Period, all service to the equipment must be provided by a Guldmann Certified Technician. A Guldmann Certified Technician is a technician who has successfully completed Guldmann Service Training, and who holds a valid Service Training Certificate from Guldmann, and is in possession of a valid password to access Guldmann's Service and Information Console (SIC). A Guldmann Service Training Certificate and SIC password are valid for three years (only USA) from the date the technician is first certified. Thereafter, the technician must undergo re-certification training to obtain a new valid certificate and password. Any parts or components repaired or replaced by a Guldmann Certified Technician will be guaranteed for the remainder of the Warranty Period. In the event the warranty is rendered null and void, the purchaser shall indemnify and hold Guldmann harmless of and from any and all claims or liability arising as a result of equipment malfunction or misuse.

---

#### **B. Service or Repair**

Contact Guldmann Repair for an authorization to return any defective item during the Warranty Period. You will be provided with a return authorization number and address for returning the item for warranty service or replacement. Do not return items to Guldmann under warranty without receiving a Return Authorization Number.

If mailing the item, pack it carefully in a sturdy carton to prevent damage. Include your Return Authorization Number, a brief description of the problem and your return address and phone number. Guldmann does not assume the risk of loss or damage while in transit, so it is recommended you insure the package.


**8.00**

**Fitting Multi support sling – lifting upper body**

Place the sling at the upper part of the lumbar area and lower part of the upper body.


Place palm on the sling and press down on the underlay while pushing the sling in underneath the patient.

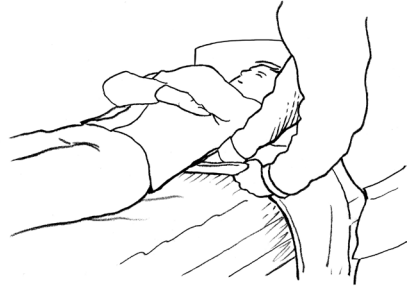
Pull the last piece of the sling out on the other side so that the sling is equally distributed on both sides of the patient.

Attach the lifting straps to the lifting hanger, and lift the patient's upper body off of the underlay.

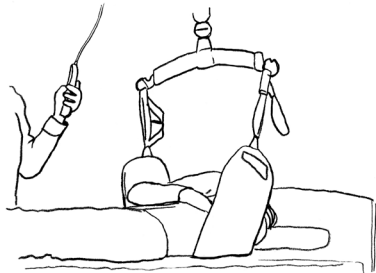
Tuck the sling in under the leg, where it is most appropriate for the care task to be performed.

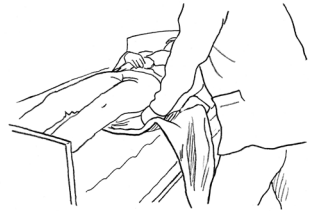
Attach the lifting straps to the lifting hanger.

To lift the patient's leg while his/her knee is bent, place the sling under the top half of the lower leg.

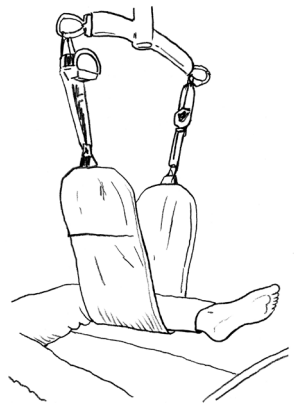To lift the patient's leg while his/her knee is extended, place the sling under the bottom half of the lower leg. If necessary, support the patient's thigh with a pillow or a bend in the bed.

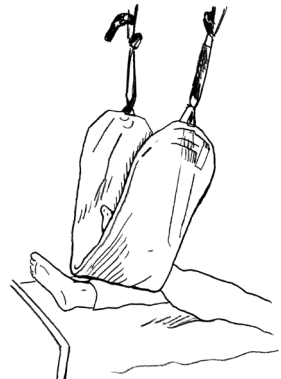**10.00**

**Fitting Multi support sling – lifting leg from wheelchair**

Lead the sling in underneath the patient's thigh and the back side of his/her knee.

Carefully push and pull the sling in underneath the patient's leg so that the sling is equally distributed on both sides of the leg.


After the lifting strap has been attached to the lifting hanger, the leg can be lifted from the underlay.


**11.00**

**Fitting Multi support sling – other uses**

Lifting of pelvis

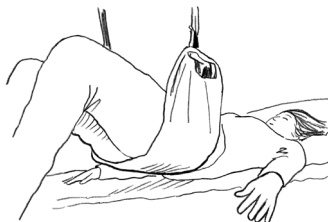
Lifting leg with patient lying on his/her side:  
(requires 2 slings)

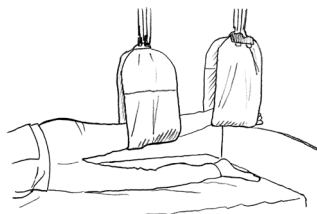---

Lifting leg of a patient who needs to get in/out of bed:

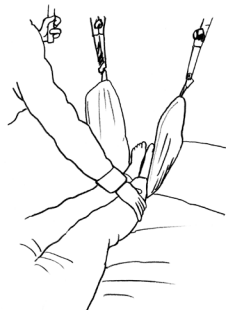
---

Lifting arm.

---

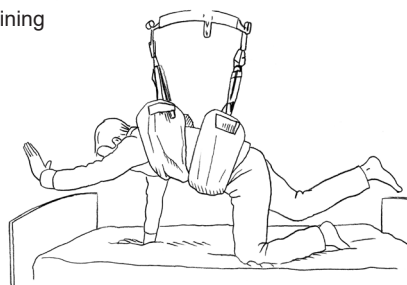Knee-four-standing balance training

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## **Guldmann GmbH**
## **Guldmann Sverige AB**
## **Guldmann Srl**
## **Guldmann Inc.**